

The Rondo logo features the word "ROND" in a bold, blue, sans-serif font, followed by a circular icon containing a stylized blue sphere with a white outline.

ROND

Cannulated screw with integral adjustable washer

- ✓ *Stable fixation*
- ✓ *Easy removal*



We invented it

citieffe[®]
Essential Moves in Trauma

Surgeons must always decide on the best approach to follow according to their clinical judgement and the patient's needs. Citieffe does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery. The information presented is intended to demonstrate the extent of Citieffe product offerings. Before use a surgeon must always consult the package insert, product label and/or instructions for use. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Citieffe representative if you have questions about the availability of Citieffe products in your area.

Cannulated screw with integral adjustable washer

Compression and adaptability to the surface of the bone


Orientable washer: the geometry of the screw head allows best fit on the bone

Versatile and intuitive system

Color coded implant and instrument set for a facilitated, rapid and safe operative technique

A wide choice of diameters (from 3mm to 8mm) and **lengths** (from 10mm to 160mm)

Partially or total threaded screws

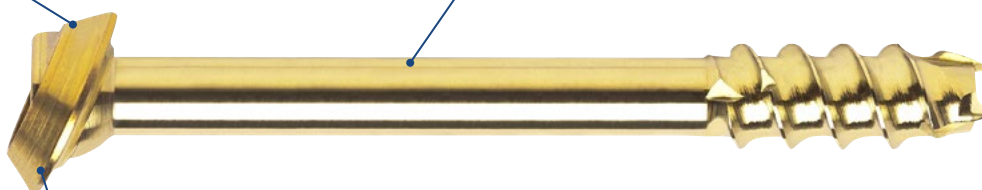


Stable fixation and easy removal


The **integrated washer** remains with the screw during insertion and removal

Single pitch thread for a greater compaction and compression of the fracture


Double pitch thread for a faster insertion speed

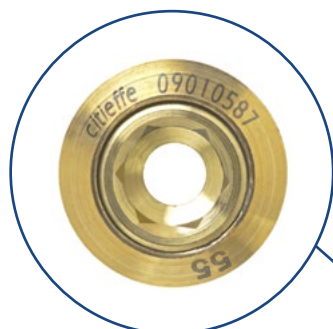

Integral washer

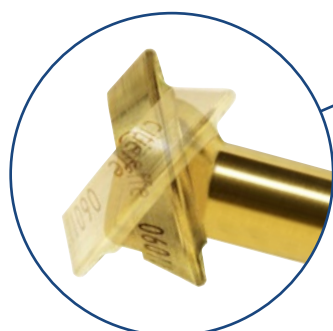

Adjustable washer

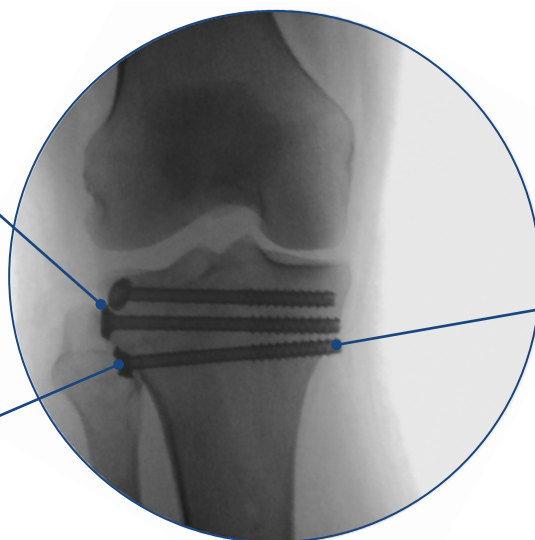

Cancellous thread

Bibliography

- Orengo G., 2008. Le viti rondò: semplicità e affidabilità di una nuova vite cannulata. *Proceedings from XIV Congresso Nazionale C.I.O.*
- Addevico F., 2019. La forza di inserzione di una vite come parametro di giudizio della sintesi. *Proceedings from 104° Congresso Nazionale S.I.O.T.*
- Addevico F., 2019. Strategie di ripristino della capacità meccanica di una vite spanata salvando la location. *Proceedings from 104° Congresso Nazionale S.I.O.T.*
- Addevico F., Morandi M., Scaglione M., Solitro G.F., 2020. Screw insertion torque as parameter to judge the fixation. Assessment of torque and pull-out strength in different bone densities and screw-pitches. *Clin. Biomech.* 72, 130-135

Rondò screws for every anatomical district

ø3mm

single pitch thread
variable length



ø4mm

single pitch thread
variable length



ø5mm

single pitch thread
variable length



ø6.5mm

double pitch thread
length 16mm



ø6.5mm

double pitch thread
length 22mm



ø6.5mm

double pitch thread
length 32mm



ø7mm

single pitch thread
length 20mm



ø7mm

single pitch thread
length 40mm



ø7mm

single pitch thread
total threaded



ø7.5mm

double pitch thread
length 20mm



ø7.5mm

double pitch thread
length 35mm



ø8mm

single pitch thread
length 20mm




ø8mm

single pitch thread
length 40mm



ø8mm

single pitch thread
total threaded

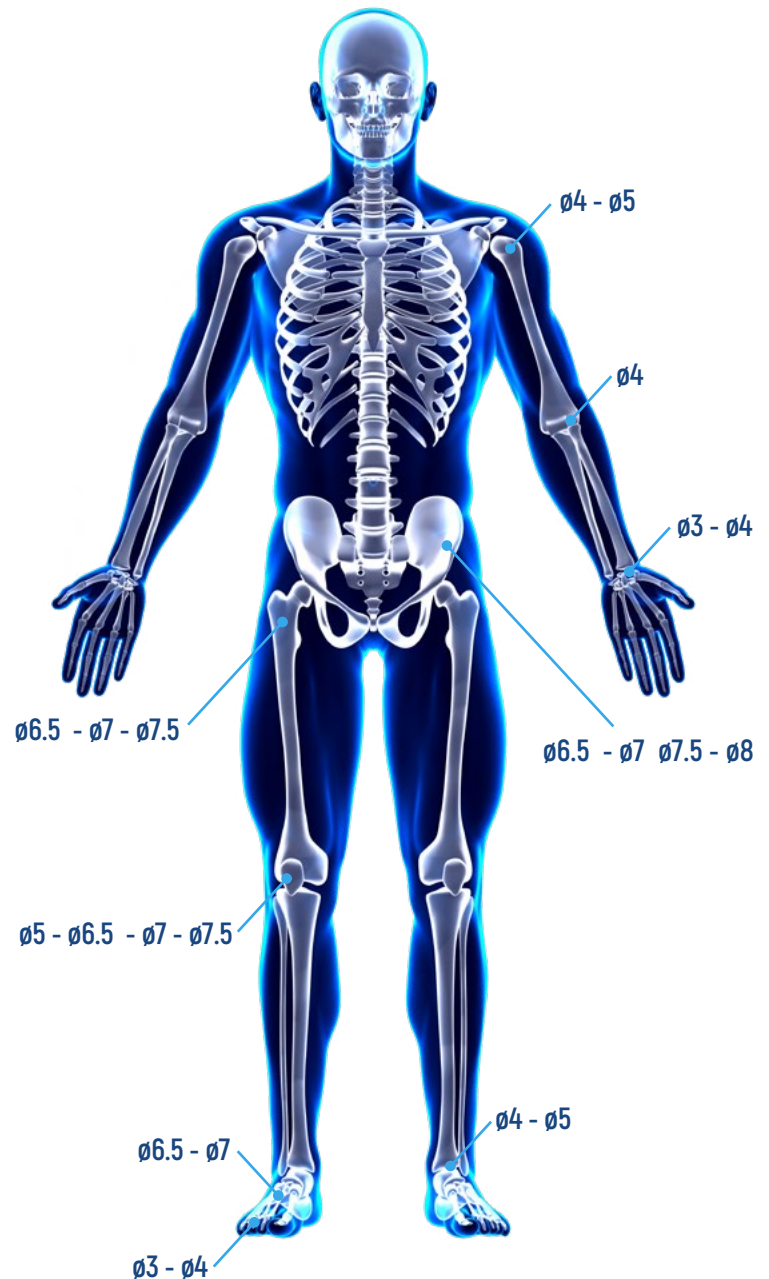



Cannulated screw with integral orientable washer



**INSPIRED BY PEOPLE
MOVED BY INNOVATION**



Citieffe s.r.l.

Via Armaroli, 21

40012 Calderara di Reno (BO) - Italy

Tel +39 051 721850 - Fax +39 051 721870

info@citieffe.com - www.citieffe.com

YOUR DISTRIBUTOR IS