


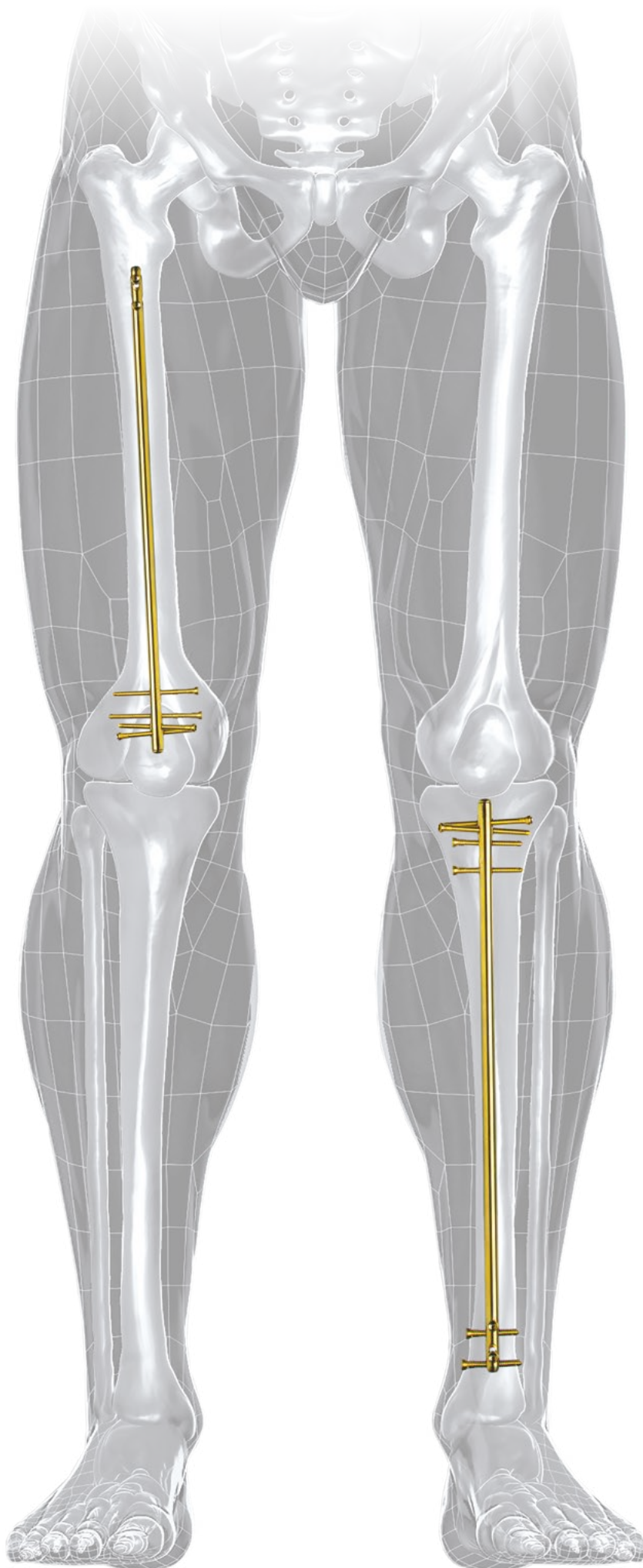

ESTREMO
CITIEFFE NAILING SYSTEM

Lower limb nailing system

✓ Unique design
for lower limb

✓ Multiple
indications


✓ Versatile
locking options

citieffe[®]
Essential moves in Trauma

unique design

Simplified suprapatellar approach

Operative technique performed entirely through the soft-tissue cannula

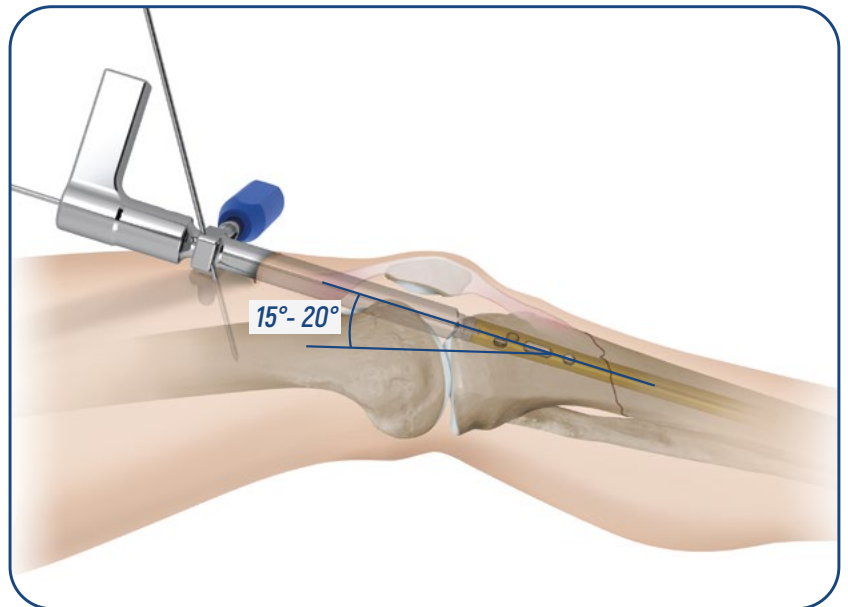


- ✓ Anatomical bend for an easy suprapatellar nail insertion
- ✓ Different multidirectional screws locking options
- ✓ Multiple static and dynamic locking options

Multi-directional proximal locking

- 2 oblique screws (40° angle)
- 2 parallel screws with mediolateral approach

Tibial plateau bone fragment reduction

- 2 guided out-nail screws

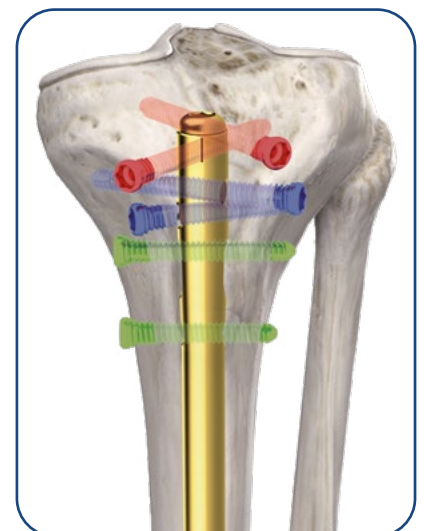

designing multiple indications

slot for end cap
and nail removal

2 holes for
oblique screws

slot for dynamic locking
and intraoperative compression

hole for static locking

Nail

- ✓ Titanium alloy ASTM F136
- ✓ Solid nail: $\varnothing 8\text{mm}$
- ✓ Cannulated nail: $\varnothing 9, 10, 11, 12, 13\text{mm}$
- ✓ Available lengths: 150mm, 200mm, from 260mm to 470mm (with 15mm increments)

Anti pull-out system for the oblique screws



Screws



$\varnothing 5.2\text{mm}$ cortical screw




$\varnothing 4\text{mm}$ cortical screw (for solid nail)



$\varnothing 5.2\text{mm}$ cancellous screw

Multiple options to enhance the stability

Multidirectional static and dynamic locking of screws

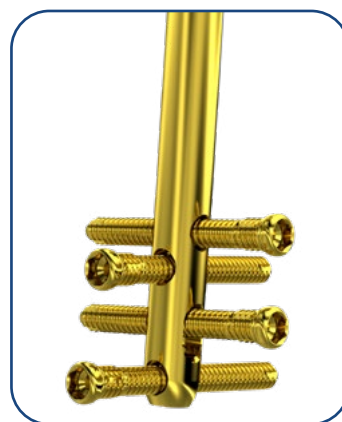


3 static locking options

slot for static or
dynamic locking


indications

Easy retrograde approach

Operative technique performed entirely through the soft-tissue cannula

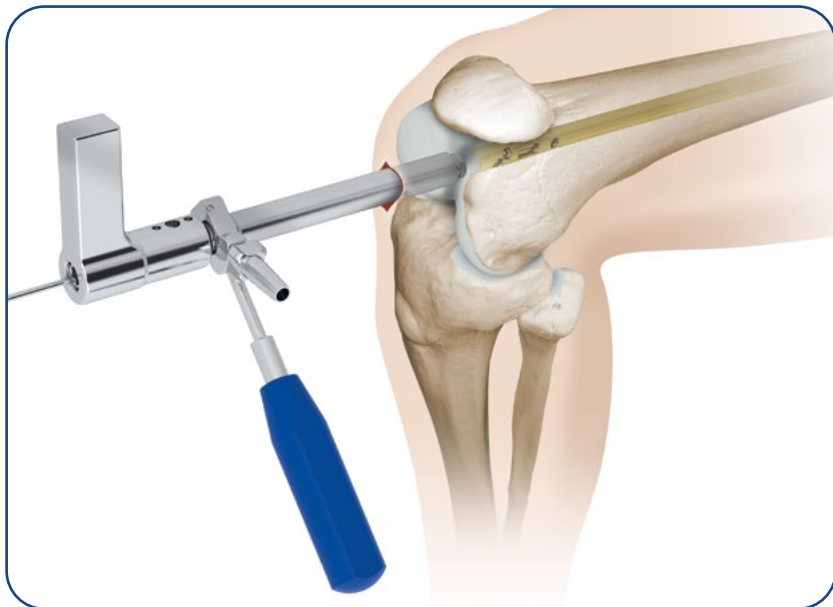

- ✓ Ideal radius of curvature for retrograde approach
- ✓ Universal nail for right and left limb
- ✓ Guided screws locking

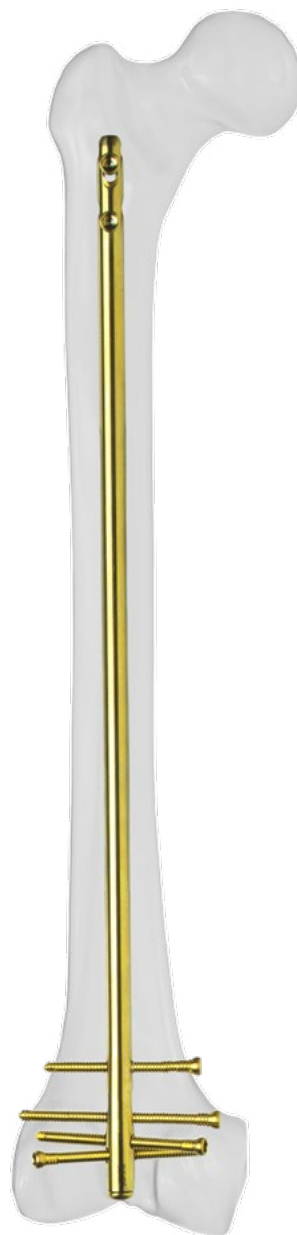

Single instrument set for all approaches



Centering guide for an easier insertion of the out-nail screws



Patented cannula locking system



Screwdriver **with retention system** for all the screws and caps

ESTREMO
CITIEFFE NAILING SYSTEM

Lower limb nailing system

citieffe® | INSPIRED BY PEOPLE
Essential moves in Trauma | MOVED BY INNOVATION



Citieffe s.r.l.

Via Armaroli, 21

40012 Calderara di Reno (BO) - Italy

Tel +39 051 721850 - Fax +39 051 721870

info@citieffe.com - www.citieffe.com

YOUR DISTRIBUTOR IS